


visicolor® ECO Ammonium 15


Ammonium 15/Ammonium 15/Amonio 15

Method(e) / Método		PF-12 585 nm	PF-12 ^{Plus} 585 nm	Advance 585 nm
5101	0.4–6.2 mg/L NH ₄ -N	■	■	■
5102	0.5–8.0 mg/L NH ₄ ⁺	■	■	■
5103	0.5–8.0 mg/L NH ₃	■	■	■

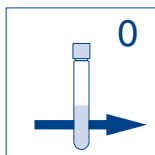


Null mit Probe
Blank with sample
Blanc avec échantillon
En blanco con muestra


Geeignet für Meerwasser nach Verdünnung 1+9
Suitable for sea water after dilution 1+9
Convient à l'eau de mer après dilution 1+9
Adecuado para el agua de mar tras su dilución 1+9

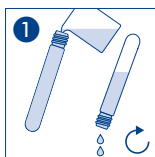


Küvetten spülen
Rinse test tube
Rincer les éprouvettes
Lavar los tubos

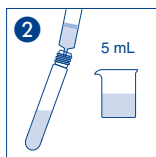


Probe
Sample
Echantillon
Muestra

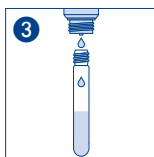


10 ∅ NH₄-1

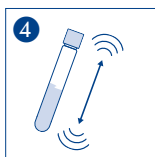

Schütteln
Shake
Agiter
Agitar

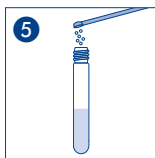


1 ∅ NH₄-2

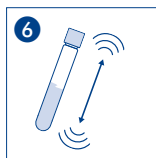


Schütteln
Shake
Agiter
Agitar




Warten
Wait
Attendre
Espere

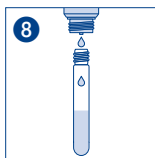

4 ∅ NH₄-3

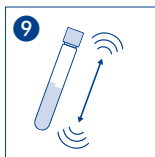


Schütteln
Shake
Agiter
Agitar

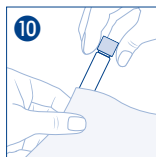

Säubern
Clean
Nettoyer
Limpiar



Warten
Wait
Attendre
Espere

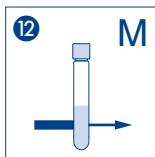

Messung
Measurement
Mesure
Medición