


visicolor[®] ECO Chrom(VI)


Chromium(VI) / Chrome(VI) / Cromo(VI)

Method(e) / Método	PF-12 540 nm	PF-12 ^{Plus} 540 nm	Advance 540 nm
5201 0.02–0.50 mg/L Cr(VI)	■	■	■
5202 0.04–1.00 mg/L CrO ₄ ²⁻	■	■	■

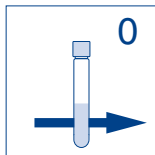


Null mit Probe
Blank with sample
Blanc avec échantillon
En blanco con muestra


Geeignet für Meerwasser
Suitable for seawater
Convient à l'eau de mer
Adecuado para el agua de mar

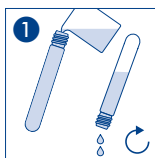


Küvetten spülen
Rinse test tube
Rincer les éprouvettes
Lavar los tubos

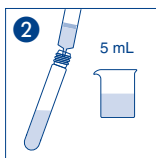


Probe
Sample
Echantillon
Muestra

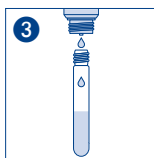


5 Δ Cr-1

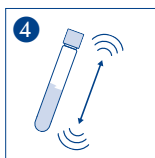

Schütteln
Shake
Agiter
Agitar

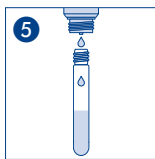

5 Δ Cr-2

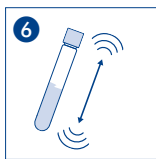


Schütteln
Shake
Agiter
Agitar



Säubern
Clean
Nettoyer
Limpiar



Warten
Wait
Attendre
Espere

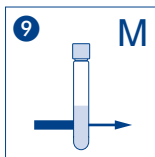

Messung
Measurement
Mesure
Medición