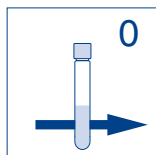


visicolor® ECO Zink


Zinc / Zinc / Zinc

Method(e) / Método


5981 0.1–3.0 mg/L Zn²⁺PF-12
620 nmPF-12^{Plus}
620 nmAdvance
620 nm

Null mit Probe
Blank with sample
Blanc avec échantillon
En blanco con muestra


Geeignet für Meerwasser nach Verdünnung 1+9
Suitable for sea water after dilution 1+9
Convient à l'eau de mer après dilution 1+9
Adecuado para el agua de mar tras su dilución 1+9

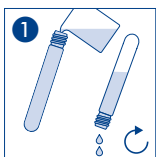


Küvetten spülen
Rinse test tube
Rincer les éprouvettes
Lavar los tubos

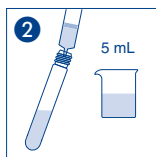


Probe
Sample
Echantillon
Muestra

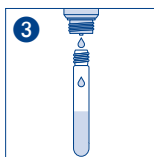


5 ♡ Zn-1

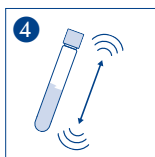


Schütteln
Shake
Agiter
Agitar

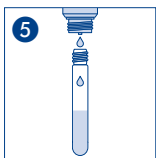


5 ♡ Zn-2

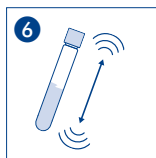

Schütteln
Shake
Agiter
Agitar

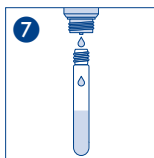


10 ♡ Zn-3

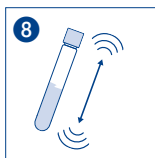


Schütteln
Shake
Agiter
Agitar



Säubern
Clean
Nettoyer
Limpiar

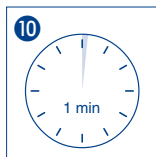

Warten
Wait
Attendre
Espere

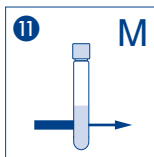

Messung
Measurement
Mesure
Medición